



SAFETY TRAY

ENZYMATIC CLEANER



Instruments







Monarch Enzymatic Cleaner is formulated with 5 active enzymes and the latest in advanced surfactant synergy, helping you to properly clean your instruments prior to sterilization.

## Monarch Enzymatic Cleaner

- For use as presoak cleaner, ultrasonic detergent and for manual cleaning.
- Rapid removal of debris resulting in labor savings.
- Disintegrates proteinaceous and fat materials, sugars and carbs in the tiniest crevices.
- Neutral pH, non-corrosive resulting in protection of valuable instruments for longer life.
- Non-foaming.
- Available in convenient unit-dose packets and 84.5 oz. bottle.
- Mixes instantly.
- Fresh scent.

Monarch Enzymatic Cleaner contains a neutral pH formula and higher concentrations of detergents and enzymes, resulting in:

- Optimal cleaning results.
- Protection of valuable instruments for longer life.
- Process efficiencies.

This uniquely formulated, presoak cleaner and ultrasonic detergent will effectively remove blood, tissue, fats, carbs, sugars, and other difficult debris from surgical instruments, and assist in complete cleaning before sterilization.



Soak and hold instruments in Monarch Safety Tray to make them easier to clean later. The tray also provides a safer way of transporting the instruments from the operatory to the instrument processing area.


# Monarch Safety Tray

- Facilitates easy and safe instruments processing.
- Removable perforated tray allows for easy rinsing, draining and pick-up of instruments.
- Hinged cover protects staff from bath-generated aerosols.
- Constructed of heavy gauge, high-impact plastic.
- Holds 3 quarts; fill lines designate 1-quart increments.

### Surfactant Synergy

Monarch Enzymatic Cleaner contains anionic and nonionic surfactants. These surfactants ensure that the instruments are clean and free of fragments that can redeposit, especially in presoaking, manual and ultrasonic cleaning, and are excellent in removing lipid/fat soils.

#### Instrument Processing Cycle

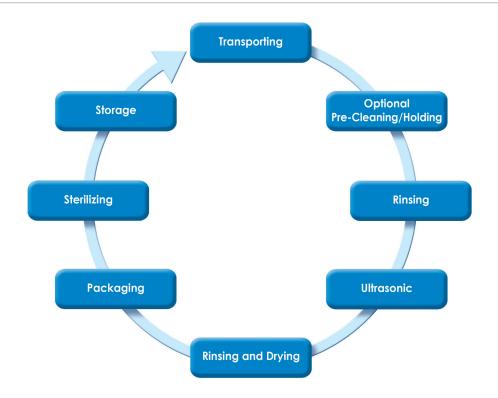

#### airtechniques.com

#### Order Information

Order No.	Product	Description
H6310-N	Monarch Enzymatic Cleaner	Box of 50 unit-dose packets; 1 packet makes 1 gallon (12 boxes per shipper case)
H6201	Monarch Enzymatic Cleaner	84.5 oz. bottle; makes 125 gallons (4 bottles per shipper case)
H6150	Monarch Safety Tray	13" L x 8.25" D x 5.1" H; holds 3 quarts





CORPORATE HEADQUARTERS 1295 Walt Whitman Road Melville, New York 11747- 3062 800-247-8324



Monarch is a registered trademark of Air Techniques, Inc. © Copyright 2018 PN 9922-29 RevC

