



CLEANSTREAM & DISPENSER

Evacuation System Cleaner & Dispenser with **VORTEX** TECHNOLOGY

Optimizes suction performance to keep your vacuum running at its best.



CleanStream
by
AIR equipped for life
TECHNIQUES

RELIABLE TEAMWORK

CleanStream is the ONLY evacuation line cleaner Air Techniques recommends using daily to ensure highest performance and longevity of your equipment.

Superior cleaning effectiveness

- A proprietary blend of surfactants and chelating agents, removes debris and prevents future build-up maintaining suction performance to ensure an uninterrupted workflow.
- Ideal cleaning results are achieved when CleanStream is used together with the Dispenser with Vortex technology.
- Works with both wet and dry vacuums.

Safe for you and your equipment

- Cleaning suction lines daily with Vortex action ensures optimal aerosol mitigation during procedures.
- Non-foaming solution protects vacuum for longer life.
- Non-corrosive & non-hazardous.
- Amalgam separator safe with a neutral pH.

Cost-Effective & Convenient

- Only 20 ml of the concentrate is needed to clean 1 operator and deodorize lines.
- Solution stays effective for up to 5 days when premixed.



Clean and protect dental office plumbing lines and suction tubing

Don't get caught off guard with an unexpected shutdown due to clogged evacuation lines in your facility. Keep them clean with a state-of-the-art system utilizing CleanStream Dispenser with Vortex Technology.



Without CleanStream & Vortex Technology



With CleanStream & Vortex Technology

Hands-Free Cleaning Made Simple

The **CleanStream Dispenser with Vortex Technology** enables hands-free evacuation-line cleaning that is easy, consistent, and mess-free—with no spills or splashes.

Protect Your Evacuation Lines

Help prevent buildup and unexpected blockages that can lead to reduced suction, equipment issues, and costly downtime. Daily cleaning keeps your evacuation lines flowing freely and your practice running smoothly.

Support Effective Aerosol Capture

Clean evacuation lines help maintain the airflow your HVE system needs to effectively capture aerosols during procedures. Make daily line cleaning part of your practice's safety routine.

Help Extend Vacuum System Life

Routine maintenance helps protect your investment. Clean your evacuation lines daily with **CleanStream and Vortex Technology** to help keep your vacuum system performing at its best.

Exact Dispensing Every Time


Precisely dispenses the exact amount of solution needed to clean one operator. No waste. No guesswork.



1st operator: Place dispenser upright. Half the solution aspirates automatically, then stops.



2nd operator: Turn dispenser on its side to aspirate the remaining solution.

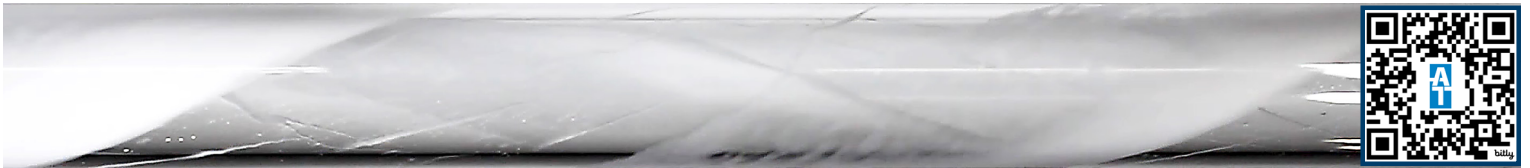

To see Vortex Technology in action, click link or scan QR code: <https://youtu.be/Z84Z5BZLgio>



Steven H. Berk, D.D.S.
Annandale, VA

“My assistants
enjoy the ease
and efficiency
of how the
system works.”

5 Years as the Preferred Product

Order Information

Order No.	Product	Description
57660	CleanStream Starter Kit	Includes: 1 Dispenser and 1 bottle (33.8 oz./1 L) of CleanStream
57850	CleanStream Evacuation System Cleaner	84.5 oz. (2.5 L) Bottle with integrated measuring cap - 125 cleanings
57840	CleanStream Evacuation System Cleaner	33.8 oz. (1 L) Bottle with integrated measuring cap – 50 cleanings
57083	CleanStream Evacuation System Cleaner	338 oz. (10 L) Container – 500 cleanings
57665	CleanStream Dispenser with Vortex Technology	1 Dispenser with adapters* for all suction tubings
H6210	Dosing Pump (20 ml)	For CleanStream 338 oz. (10 L) Container
H6209	Dosing Aid (Canister Tap)	For CleanStream 338 oz. (10 L) Container
H6214	Dosing Pump (20 ml)	For CleanStream 84.5 oz. (2.5 L) Container
57655	CleanStream Dispenser Cap	Replacement cap w/updated* adapters for CleanStream Dispenser

* adapters are now affixed for dispensers manufactured as of 1/2024



Download Free Utility Room
Maintenance Checklist



CORPORATE HEADQUARTERS
1295 Walt Whitman Road
Melville, New York 11747- 3062
800-247-8324



© Copyright 2026
PN 9922-633 Rev A

